

Soprano Titanium

Stronger, Faster, Smarter

Hair Removal, Reinvented.

Introducing Soprano Titanium

A breakthrough hair removal platform that combines ultimate functionality with unparalleled comfort. Thanks to its large 4 cm² spot size and advanced cooling system, treatments are now much faster and pain free – offering a significantly improved solution for patients and practitioners alike. Soprano Titanium leverages the benefits of three combined wavelengths, with an improved patient experience and a business-oriented approach, creating a unique and result-driven new solution in the world of professional hair removal.

The Titanium's **System Highlights**

Exceptionally Large Spot Size

40% faster!

Soprano Titanium features TRIO 4cm², a unique and efficient new applicator with an especially large spot size of 4 cm², capable of working on a grid of 600cm square. This breakthrough applicator makes treatments much faster and pain free, significantly cutting treatment times and ensuring a more comfortable experience for patients.

Dual Connector

More Flexibility. Less Fuss.

Soprano Titanium is the first version of Soprano hair removal systems that includes two connectors, which enables two different applicators to be connected to the platform at all times. This feature allows for more flexibility, while eliminating time spent on switching applicators during treatments.

ICEPlus^{*}

Continuous Contact Cooling

ICE Plus is an advanced technology that continuously cools skin and facilitates range-controlled temperature throughout the entire treatment. The ICE Plus[™] cold sapphire tip minimizes epidermal risks while maintaining heat within the dermis where hair follicles are treated. Temperature remains controlled and constant which contributes to the diode's efficiency. Patients can now enjoy a completely pain-free, cool and comfortable experience throughout the session.

15' Android Screen

Friendly and Functional

Equipped with a new large sized 15' Android screen, with pre-set parameters, practitioners can enjoy a more intuitive, user friendly work experience.






Smart Clinic

Manage Your Clinic. Maximize Your Potential.

Smart Clinic is a pioneering, cloud-based businesses development tool that offers continuous access to live data through an application and online platform. It empowers better management processes simply and seamlessly, reflecting a real time status of platforms, applicators, and productivity, letting you access data and make adjustments without losing time. Smart Clinic drives optimization and standardization, empowering decisions based on genuine data and statistics – not guesswork. It also boosts marketing endeavors by enabling notifications of launches, new products, event invitations and more.

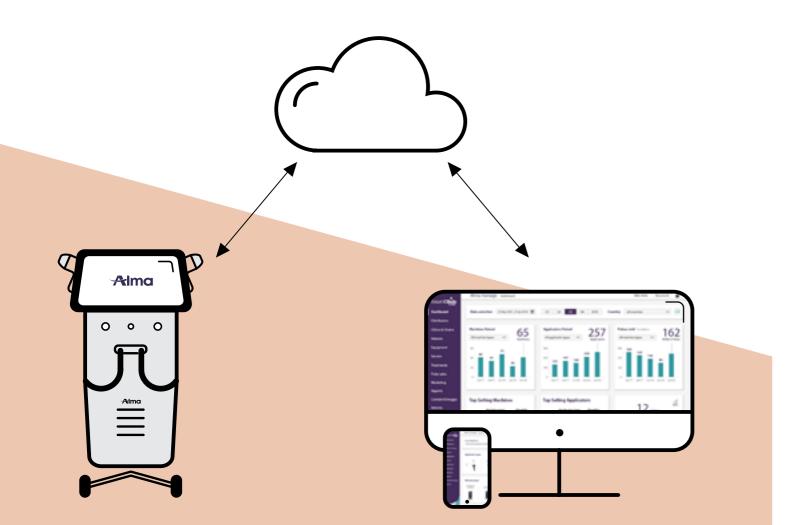

Never Miss a Beat

Smart Clinic is a multi-faceted solution that provides a complete picture of your business; whether it is a single clinic or a chain.

Platforms and Applicators | Lets you view the number of pulses, time operated, areas treated, number of platforms and applicators operated by each practitioner at any given time. This data is crucial for ongoing modifications of your clinic's equipment and for collecting operational statistics (per platforms, applicator and clinic), tracking practitioner performance and improving human resource management.

Treatment Management | Follow up on individual patient care, including areas treated, protocols, and date of visit. Monitor patient cases and needs, and contact them when their next treatment is coming up. Improve customer support and experiences, enhance clinical results and reinforce customer engagement and relationships. Smart Clinic presents clinical protocols on the screen, making hard-copy protocols a thing of the past.

Maintenance and Service | Smart Clinic lets you open up a service call at any time through the application, sending you a troubleshoot guide so that problems can be solved when they arise. It offers the ability to track treatments and view statistics for each machine to prevent overuse and underuse, follow up on errors and receive ongoing maintenance alerts (water/filters). With Smart Clinic, you can prevent downtime, manage your resources more efficiently and boost productivity.

Synergetic Technologies. Superior Results.

SHR[™] technique SHR[™] - Alma's proprietary hair removal innovation is the only clinically proven method of laser hair removal that is virtually painless, and it is the safest laser hair removal method for all skin types, including dark and tanned skin. Its unique gradual heating and rapid/fast pulse delivery method effectively damages the hair follicle and hinders re-growth, while preventing injury to the surrounding tissue. SHR[™] affects the melanin of the skin using a unique method, allowing the procedure to be performed on all skin types anytime of the year as well as on tanned skin.

3D Technology




Featuring its exclusive 3D technology, Soprano Titanium combines the three most effective laser wavelengths into a single applicator, simultaneously targeting different tissue depths and anatomical structures within the hair follicle. By combining the absorption and penetration levels of different wavelengths, along with extended treatment coverage, increased comfort and low maintenance requirements, Soprano Titanium is the safest and most efficient hair removal treatment available today. 3 dimension technology for one superior hair removal solution.

ALEX 755nm The Alexandrite wavelength improves the absorption of powerful energy by the melanin chromophore, making it ideal for the most extensive range of hair colors and types, especially light-colored and thin hair. With increased superficial penetration, the 755nm wavelength targets the Bulge of the hair follicle and is especially effective for superficially embedded hair.

SPEED 810nm The 810nm is a classic hair removal wavelength, offering deep penetration of the hair follicle with high average power, and a rapid repetition rate and a large 2cm². spot size for time-efficient treatments. The 810nm has a moderate melanin absorption level, making it safe for darker skin types. Its deep penetration capabilities target the Bulge and Bulb of the hair follicle, while moderate tissue depth penetration makes it ideal for treating the arms, legs, cheeks and chin.

YAG 1064nm The YAG 1064 wavelength is characterized by lower melanin absorption, making it a perfect solution for darker skin types. The1064nm also offers the deepest penetration of the hair follicle, so that it targets the Bulb and Papilla, while deeply treating embedded hair in areas such as the scalp, armpits and pubic areas. With higher water absorption generating a higher temperature, incorporating the 1064nm wavelength increases the thermal profile of the overall treatment, leading to highly effective hair removal and improved results.

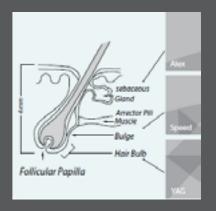

The Titanium's **Applicators**



Complete more treatment sessions in less time with this large spot size applicator. Treat larger areas faster, while offering patients a safe and highly effective hair removal procedure.



TRICA.

3D and 40% faster!

This exclusive applicator is a valuable asset for any clinic. With its large 4 cm² spot size, and integration of three wavelengths, it facilitates increased treatment coverage to significantly reduce treatment time and make the process easier and safer for patients. With TRIO 4cm², hair removal is quick, efficient and convenient.



COMPACT 810 - Light And Lean

With its ergonomic, light weight design, this applicator makes treatments much more comfortable for practitioners and more enjoyable for patients.



TRIO - Three In One

This applicator enables the delivery of three highly effective wavelengths in one. Covering a spot size of 2cm², it simultaneously targets different tissue depths and anatomical structures within the hair follicle, increasing treatment efficiency.



DIODE LASER ALEX 755 -The Power Of Diversity

Powerful absorption, combined with continuous cooling, treatment coverage, comfort and low maintenance, enables the treatment of a wide range of hair types and colors, especially light-colored and thin hair.



FACIAL TIP - Improved Reach

Easily treat areas that are normally hard to reach, including the ears, nostrils and glabella.

Platform $\left(\begin{array}{c} \hline \\ \hline \\ \hline \\ \hline \\ \hline \\ \end{array} \right)$ **Benefits**

*

3D

(모)

I-IV

LARGE 4 cm² SPOT SIZE FOR QUICK TREATMENTS ICE PLUS[™] CONTINUOUS COOLING SYSTEM TWO CONNECTORS SMART CLINIC DATA CENTER

CLINICALLY PROVEN 3D TECHNOLOGY

INTUITIVE, USER FRIENDLY LARGE SCREEN (15' ANDROID) WITH PRE-SET PARAMETERS

I-VI SKIN TONES AND TANNED SKIN

NO CONSUMABLES

I am very pleased with the results I am able to achieve with the Soprano Titanium. Besides for the high efficiency and safety record, which are the base for any hair removal treatment, there are additional features, such as Smart Clinic and the unique TRIO 4cm² applicator, which set it above and beyond any other system on the market today.

Dr. Pablo Naranjo MD, PhD, Medical Director of Elite Laser Clinic & Laser Unit at NISA Hospital, Madrid, Spain



Soprano Titanium Hair Removal, Reinvented.



Alma Lasers GmbH

Nordostpark 100-102 90411 Nuremberg, Germany Tel. + 49 911 / 89 11 29-0 Fax + 49 911 / 89 11 29-99

Email: info@almalasers.com

www.almalasers.com

© 2020 Alma.All rights reserved. Alma ,its logo, Soprano Titanium, TRIO, ICE Plus, Smart Clinic, SHR & Speed are trademarks and/or registered trademarks of Alma in the United States and/or other countries.

PBSP28111801_03



